

Non-human Primate Startup, Shutdown and Decontamination Protocol

Please note: HTS is **NOT** permitted for NHP work.

Startup


1. Empty the waste tank at the start of every NHP run and then fill the waste tank with 1500mL of bleach to ensure a minimum concentration of 10% after your run– this will ensure you know how much bleach is present at the start of your run.
2. Continue with normal startup for Fortessa

Shutdown

1. Place a tube of 3mL 15% bleach with the loader arm to the side for 1.5mins (this should all but empty the tube – if it does not, find JVD as there could be a clog in the droplet containment lines that need to be addressed.
2. Run on “High” with a new 15% bleach tube for 15 mins.
3. Place a tube of 3mL 10% Contrad bleach with the loader arm to the side for 1.5mins (this should all but empty the tube) – if it does not, find JVD as there could be a clog in the droplet containment lines that need to be addressed.
4. Run on “High” with a new 10% Contrad tube for 5 mins.
5. Place a tube of 3mL DI H₂O with the loader arm to the side for 1.5mins (this should all but empty the tube) – if it does not, find JVD as there could be a clog in the droplet containment lines that need to be addressed.
6. Run on “High” with a new DI Water tube for 5 mins.
7. Discard ALL cleaning tubes into the red biohazard bin.
8. Empty the bleach treated waste from the waste tank
9. Add 2 inches of bleach to the bottom of the waste tank.

Decontamination of surfaces

The user will spray Cavicide (obtained from JVD prior to run) on the surfaces indicated in the image below and allow to sit on the indicated surfaces for 5mins. After sitting wipe the surfaces with a paper towel to remove Cavicide.

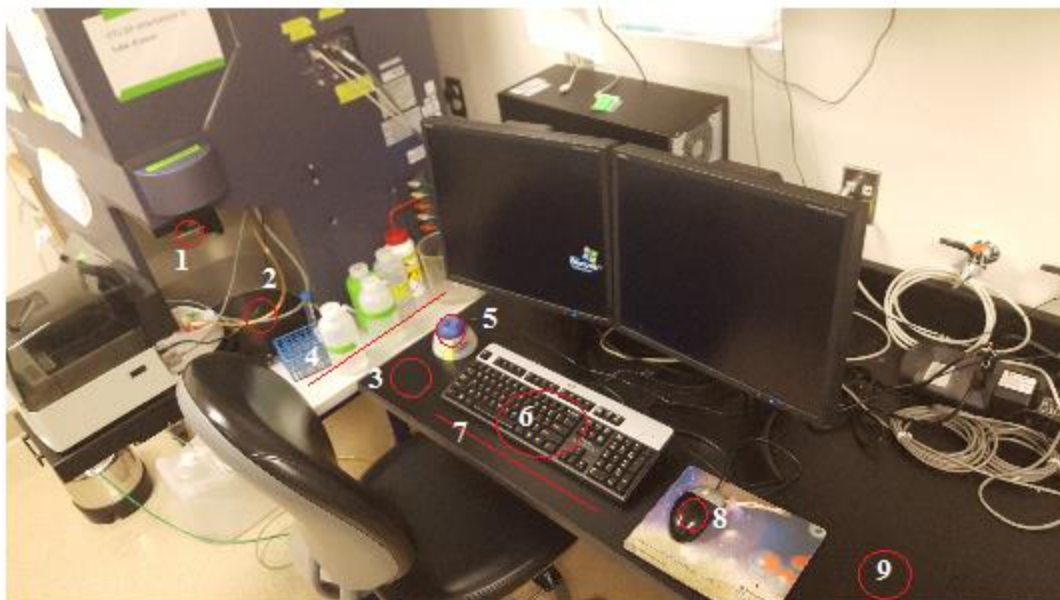

9 total surfaces

1. Sample Injection Tube
2. Fortessa stage area
3. Desk left side
4. Cleaning solution and tube rack area
5. Vortex (If used)
6. Keyboard
7. Desk- under keyboard
8. Mouse
9. Desk right side